


The Band Sling has a wide body component for increased support of the upper torso; minimizing the sling "riding up" under the arms of the person. Designed for use with BHM Medical sit and stand floor lifts, this sling is to be used for seat to seat transfers, toileting and pericare.

The sling allows for removal of lower extremity clothes of the person without removal of the sling.

It can also be used for practice weight bearing and balance for partially dependant people.

Persons should have good upper body control to hold themselves up and lean back.

The person's arms are positioned outside the sling at all times.

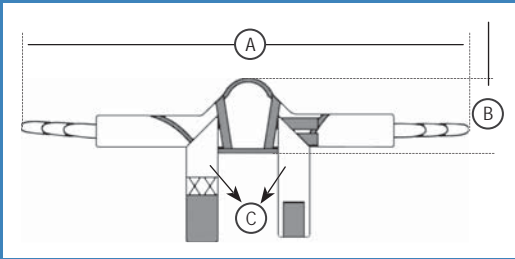

# Band Sling

# PRODUCT	DESCRIPTION
#TST-S	Small Padded Band Sling
#TST-M	Medium Padded Band Sling
#TST-L	Large Padded Band Sling
#TST-XL	Extra Large Padded Band Sling

FEATURES	
Warranty against defect	1 year - washable models*
Safe working load	272 kg / 600 lb / 42.86 st
Open area from above waist to thighs	Excellent for toileting. Allows easy access to undergarments, incontinence pads, or for cleaning. Sling supports the upper body
Padded and quilted polyester	Special designed fabric is ultra-soft and long lasting. Sling provides extra comfort
Velcro and double belted waist band	Velcro closure adjusts to different sizes. Double belted waist band secures the person in the sling. Enables people to stand with the security of knowing they are fully supported. Easy to install on a seated person
Machine washable	Easy to clean and care for
Extra padding under arms	Provides additional cushion. Flexible and gentle
Mesh Back Insert	Sure-grip material keeps the sling in place. Soft material is gentle to skin
Back support strips	Distributes the weight evenly throughout the sling. Covers a wide area of the back for support. Reduces pressure under the arms
Strong nylon/polyester straps with positioning loops	Adjustable straps allows people with large midsection. Provides flexibility for many types of people

## TECHNICAL SPECIFICATIONS

Conforms to ISO 10535 & Marked CE

\*Conditional to normal use and maintenance according to appropriate labels, the sling must be used only with the appropriate BHM Medical Inc. lifting devices.

BAND SLING				
#PRODUCT	DESCRIPTION	DIMENSIONS		
			CM	IN
#TST-S	Small Padded Band Sling Capacity 20-45 kg / 45-100 lb / 3.15-7.14 st	A	150	59
		B	31	12
		C	109	43
#TST-M	Medium Padded Band Sling Capacity 45-90 kg / 100-200 lb / 7.14-14.29 st	A	150	59
		B	33	13
		C	122	48
#TST-L	Large Padded Band Sling Capacity 90-180 kg / 200-400 lb / 14.29-28.57 st	A	158	62
		B	37	14
		C	145	57
#TST-XL	Extra Large Padded Band Sling Capacity 180-272 kg / 400-600 lb / 28.57-42.86 st	A	163	64
		B	38	15
		C	168	66

Refer to Sling Operating & Product Care Instructions Guide for complete application details.

